Invacare[®] | SB[®] 755






When comfort comes first

Combining comfort, ease of handling and safety, the *Invacare* **SB 755** is ideal for domestic and community care environments.

The **SB 755** offers optimal comfort with an ideal seating positioning to minimise shear and friction. Carers can adjust the bed to meet individual user requirements.

Thanks to a unique integrated feature, the **SB 755** can be adapted to taller users as both the head and foot of the bed can be extended. All these features provide real benefits to users and carers in terms of comfort and general well-being.

Now available with a Light Oak colour option*

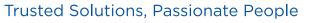







Safety is key

As always at Invacare, safety is our key concern. Side rails in both the standard bed length and the extendable version meet the increased strength demands of

Your comfort is our focus

The **SB 755** comes with several features to ensure an excellent seating position and to reduce shear and friction during movement. These include improved auto-regression, optimised length of mattress platform sections and a leg section that can be lowered below horizontal. Invacare strives to cater for all patients and so offers a length extension of 100 mm at the head and foot of this model to meet the needs of taller users.

Ease of handling

No part of the **SB 755** weighs more than 25 kg and dismountable scissors are available as an option.

Mounting and dismantling the bed has been made particularly easy to allow for greater ease of transport and storage. The bed can be stored on its base and rolled directly from van to bedroom. Assembly and transportation can be carried out by one person quickly and with minimum effort.

the IEC 60601-2-52 standard for beds, without compromising the ergonomics or working conditions of the carers.

Features and Options



Auto regression - improved functionality. Increased comfort when activating backrest, minimising shear/ friction forces.



Optimised dimensions of mattress platform according to anthropometrical measurement, for increased comfort.



Extendable mattress platform. Additional 100 mm at head and foot. Bed easily extends with 50, 100, 150 or 200 mm to suit user needs.

Features and Options



Line Side rail

Full-length side rail meets the strength and safety demand of the IEC 60601-2-52 bed standard. The Line side rail is lightweight and easy to use.



Line Extendable side rail

Full-length side rail, which can be adjusted incrementally from 2000-2200 mm.



Removable panel

On the Victoria bed end, a removable panel allows easy access from the head and foot of the bed even with the full-length side rails mounted.



SB 755 Detachable



Bed ends can be dismounted without tools.

Aniverna in the second second

The seat section opens and the two mattress support parts are divided. The mattress platform is detached from the base with a quick release mechanism.

Need to move the bed from one place to another? No lift? Narrow staircase? ...No problem!

The **SB 755 Detachable** is designed to overcome all these challenges.

This bed is so easy to mount and dismount that it can be carried out by one person; our unique dismountable system is designed for minimum effort - no need to worry about losing screws or parts!

Designed with parts that weigh no more than 25kg each, the **SB 755 Detachable** is a novel design that is unique in the market. Such easy handling ensures that the bed is easy to carry – to whatever floor.



Mattress platform is effectively released and unlocked in the seat section.



The scissor is released from the base with the quick release handle.



Activate the quick release and the scissor is detached. The handle makes it very easy to dismount.





Bed parts, mattress and accessories are stored directly on the base with special hooks; storage and transport are facilitated and the bed takes up less space.

Invacare[®] | **SB[®] 755**





Victoria bed end Detachable panel _ Available in beech. Compatible with full length side rails.



Susanne bed end Available in beech. Compatible with full length side rails.



Anita bed end

Available beech. in Compatible with full length side rails.



Cover for wooden side rail Padded cover for wooden side rails. Cover remains in place when side rail is lowered.



*Olivia bed ends

Can be used in combination with wooden side panel. Available in light oak only.



Swivel support handle The swivel support handle provides good support for the client getting to a standing position.




*Wooden side panel

be Can used in combination with Scala Basic Plus 2 side rails. Available in light oak only.



Several hand controls: Invacare Soft Control or standard with "green" button function.



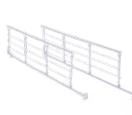



The handle can be adjusted both in height



Invacare Scala Basic Plus 2

Foldable steel side rail with plastic inserts at the mounting brackets.*



Lowering kit Using castors of 75 mm creates a lower height of 250 mm.

Mattress support

aluminium slats.

Invacare Scala Basic 2 Mattress support with Foldable steel side rail.*



Invacare Scala Medium 2 Foldable steel side rail for higher mattresses.*

Invacare Scala Decubi 2 Foldable steel side rail for extra high mattresses.*



Central brake Available as standard.

*Compliant with the new standard for beds IEC 60601-2-52.

Please see Invacare Bed Accessories brochure for more accessories.









SB755 Wide -Care solutions for plus sized people

Comprising all the features of the *SB755*, the *SB755 Wide* is available in two sizes with a mattress platform width of 105cm and 120cm. The *SB755 Wide* supports plus sized clients who require additional space for comfort and positioning. The 105cm version is particularly beneficial for carers as it reduces the risk of back injury caused by over stretching across the bed. The *SB755 Wide* has a weight limit of 200kg.

For more information, please see SB755 Wide brochure

Technical data

SB 755



950 / 1000 mm



2130 mm



adjustable:



Mattress support

dimensions:

.....



Under bed

clearance:



Backrest angle



0° - 28° Thigh angle

	Width inside: 850 / 900 mm	Length inside: 2080 mm Mattress: 2000 mm	350-820 mm Total height: 1600-2070 mm	800 - 260 - 260 - 680 mm	170 mm		gir angio
SB 755 Wide 1050 mm	Width outside: 1150mm Width inside: 1050mm	Length outside: 2130mm Length inside: 2080mm Mattress: 2000mm	Height adjustable: 350 - 820mm Total height: 1600 - 2070mm	Mattress support dimensions: 800 - 260 - 260 - 680 mm	: Under bed clearance: 150 mm	0° - 73° Backrest angle	0° - 28° Thigh angle
SB 755 Wide 1200 mm	Width outside: 1300 mm Width inside: 1200 mm	Length outside: 2130mm Length inside: 2080mm Mattress: 2000mm	Height adjustable: 350 - 820mm Total height: 1600 - 2070mm	Mattress support dimensions: 800 - 260 - 260 - 680 mm	: Under bed clearance: 170 mm	0° - 73° Backrest angle	0° - 28° Thigh angle
	, X,					Invacare Limitec	٩
	0° - 16°	0° - 18°	25° (Anale	SWL 235 kg	Heaviest part of the	Doncood Tochno	

Minimum angle Full legrest lift below horizontal

25° (Angle below frame) 14° (Angle above frame)

3 year warranty

Complies with IEC 60601-2-52

year old or smaller.

NB/ The bed must not be used by patients

SWL 235 kg Heavie (37 stone) pro Max. user weight 200 kg (31.5 stone)

under 12 years of age, or by patients with

body size equivalent to an average 12

For mattress compatibility, please refer to the user

manual which can be found at www.invacare.co.uk

Heaviest part of the product: 25 kg

Invacare Limited Pencoed Technology Park -Pencoed - Bridgend CF35 5AQ -United Kingdom Tel: +44 1 656 776 222 Fax: +44 1 656 776 220 E-mail: uk@invacare.com Sales Order E-mail: ordersuk@invacare.com www.invacare.co.uk

Invacare Ireland Ltd Unit 5 - Seatown Business Campus Seatown Road - SWORDS - County Dublin - Ireland Tel. +353 18107084 Fax +353 18107085 Email: Ireland@invacare.com www.invacare.ie

© 2016 Invacare International Sàrl All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

SB755 - UK - 04/2016 - 1523060





Washable in

a washing machine



Bariatric use

The Invacare SB755 is part of the Safe Patient Handling range.





safe**patient**handling