

HEADQUARTERS

Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T: +886 2 2268 5568
F: +886 2 2268 6525
marketing@apexmedicalcorp.com
www.apexmedicalcorp.com

PAN-ASIA

TAIWAN

Apex Medical Corp.
No. 9, Min Sheng St., Tu-Cheng,
New Taipei City, 23679, Taiwan
T: +886 2 2268 5600 / 0800 054 133
F: +886 2 2268 6525
customer@apexmedicalcorp.com

STURDY INDUSTRIAL

18, Sec.1/ Zhongxing Road, Wugu District,
New Taipei City 24872, Taiwan
T: +886 2 2981 5245
F: +886 2 2084 8208
foreign@sturdy.com.tw
www.sturdy.com.tw

CHINA

Apex Medical (Kunshan) Corp.
No. 1368, Zizhu Road, Kunshan,
Jiang Su Province, China
T: +86 512 8617 1660 Ext1150
F: +86 512 8617 8663
customer@apexmedical.com.cn

INDIA

Apex Medicalcorp India Pvt. Ltd.
303, Third Floor Block 4A DLF
Corporate Park Phase - ill,
Gurgaon, Haryana
T: +91 11 6654 4138
F: +91 11 6654 4052
enquiries.india@apexmedicalcorp.com

THAILAND

Apex Medical Thailand
65/100 Chamnanphenjati Business Center
Bldg.
11/F Rama 9 Road, Huaykwang
Bangkok 10320
T: +662 948 4848
F: +662 948 4170
sales.thai@apexmedicalcorp.com

MIDDLE EAST • AFRICA

29, Misr Liltameer Building, Stree#10, AREA#8,
Sheraton Building, Heliopolis Cairo, Egypt.
T: +202 2 2671303
F: +202 2 2678649
raafat.wahby@apexmedicalcorp.com

PAN-EUROPE

SPAIN

Apex Medical S.L.
Elcano 9, 6a planta 48008
Bilbao, Vizcaya, Spain
T: +34 94 4706408
F: +34 94 4706409
info@apexmedical.es

UK

Apex Medical Limited
Unit 33, Great Western Business Park
McKenzie Way, Worcester, WR4 9GN
T: +44 (0)1905 774 695
sales@apexmedicalcorp.co.uk
www.apexmedicalcorp.co.uk

FRANCE

Apex Medical France SAS
8 Blvd. de l'Épervière, 49000 Écouflant,
Pays De La Loire, France
T: +33 789 229 998 / +33 786 888 135
info.france@apexmedicalcorp.com

GERMANY

SLK Vertriebsgesellschaft mbH
Lindenhorster Straße 38-40 44147
Dortmund
T: +49 (0)231-925360-175
F: +49 (0)231-925360-29
info@slk-gmbh.de
www.slk-gmbh.de

AMERICAS

USA

Apex Medical USA Corp.
615 N. Berry St. Suite D, Brea, CA 92821, USA
T: +1 714 671 3818
F: +1 714 671 9040
sales.usa@apexmedicalcorp.com



Rota-Pro®

Rotational chair bed

Distributed by



Rota-Pro®

Rotational chair bed

The Rota-Pro® rotational chair bed helps people get into and out of bed without assistance, helping them to restore their mobility and independence.

Compared to alternative solutions such as hoists or care assistance, Rota-Pro® offers a simple, comfortable, and dignified means of getting into and out of bed. It comes complete with a pressure reducing mattress and bespoke bedding to hold a pillow in place. It can be configured to rotate out to the left or the right hand side according to user preference, and all positions are re-programmable for maximum comfort and support.



Included as standard

- Battery backup
- Profiled pressure reducing mattress and foot block, complete with waterproof, 4-way stretch vapour-permeable cover for comfort and ease of cleaning
- 2 x Bespoke bed sheets for mattress and foot block, including pillow retainer
- Simple, clearly labelled, luminous user hand control
- Service hand control for re-programming of profiling, rotation and seating positions to user preference
- 2 x powder coated hand rails
- Lifting pole and grab handle
- Delivery and installation (for direct customers only)
- Safety adaptor

Features/technical overview:

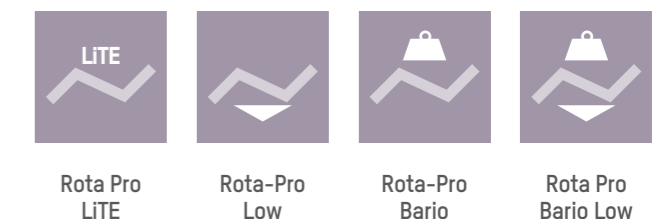
- Bed-to-chair function allows unassisted transfer to and from bed
- Chair rotation configurable to the left or the right
- Also profiles like a conventional profiling bed
- Leg-lift adjustable from nominal 0° to 40°
- Electrical adjustment of knee-break from horizontal to +40° & from horizontal to -35°
- Backrest adjustable from nominal 0° to 80° and headrest by a further 30°
- 5-section slatted mattress support
- Powder coated steel construction for durability and ease of cleaning
- Simple lift-and-lock timber side-rails
- Head and footboard and side-rails in attractive beech timber finish
- 4 x individually-braked twin wheeled castors
- Battery backup included as standard
- Unique serial number on every bed for traceability
- Linear actuators ingress protected to IPX4; Hand Controls to IPX4 for high-level protection against fluid and dust ingress
- CE marked to current EU standards
- 3 year warranty
- One button chair position
- Scissor lift frame
- Acoustic power level <50 dB(A)

Technical Parameters

Electrical specifications:	
Operating voltage	: 230V / 50 / 60 Hz
IEC 601 - 1 Device type B	
Safety class II	
Hand control	: IPX4
Main control unit	: IPX6
Power unit	: IPX4
Actuators	: IPX4

Technical Specifications:	
Safe working load (SWL)	: 170 KG 27 Stone
Maximum user weight	: 138 KG 22 Stone
Maximum height	: 83 cm
Low height	: 40 cm
Lying surface width	: 90 cm
Lying surface length	: 210 cm
Leg rest adjustment angle	: 0 to 40 degrees
Back-rest adjustment angle	: 0 to 80 degrees
Bed weight	: 192 KG

Also available in the Rota-Pro Range

Rota Pro
LITERota-Pro
LowRota-Pro
BarioRota Pro
Bario LowThree quarter length
safety hand rails

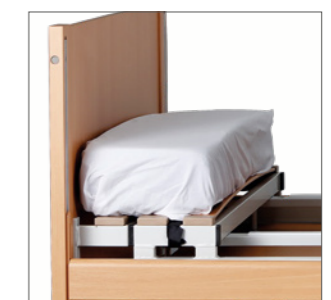
Grab handle

Full length wooden
safety sides

Intuitive hand set



Service hand control

Fitted sheets with
integrated pillow caseIndividually-braked twin
wheeled castorsAdjustable footboard with
footblock